

Protective Visor by 3DVerkstan

Published on Mar 26, 2020, By erikcederb.

License: Commercial use is allowed, you must attribute the creator, you may remix this work and the remixed work should be made available under this license.



http://www.youmagine.com/designs/26808


Hosted on





PACKAGE CONTENTS

- **1**. 6-hole clear plastic shield US letter size (8.5" x 11")
- 2. Plastic head band
- Instruction sheet and 6-hole shield template printed in the back of this instruction sheet

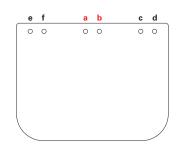



ASSEMBLY INSTRUCTIONS

Attach the 6-hole clear plastic shield (1) to the plastic head band (2): first push the two center holes **a** and **b** of the clear plastic shield into buttons **a** and **b** of the plastic head band. Once **a** and **b** holes are in place wrap the clear plastic shield carefully around the frame of the plastic head band and insert into side buttons **c** and **d** being careful not to rip the plastic when putting it on button **d** which has a slight hook feature to hold it in place. Then repeat for side buttons **e** and **f**.

TO MAKE ADDITIONAL 6-HOLE CLEAR PLASTIC SHIELDS (8.5" X 11")

MATERIALS

- 1. 3 Hole "Letter" size punch.
- 2. US Letter (8.5" x 11") clear plastic acetate (~0.2mm PET recommended).
- 3. Template: see back of this instruction sheet for template or go to https://www.youmagine.com/de signs/protective-visor-by-3dverkstan to download PDF: Laser-cut Shield North America 6 Hole Letter Format.
- 4. Scissors.

HOW TO DO THE PUNCHING

- 1. Set your hole punch to the setting for a Letter sheet.
- 2. Set or make a 1/2" spacer to create a 6-hole pattern with a 1" offset on your hole punch. The 1/2" spacer is what creates the hole offset needed. Figure A shows 1/2" spacer shown in red.
- 3. Make one punch along one of the long sides, with the short edge resting and aligned with your spacer piece.
- **4.** Flip the clear plastic sheet, and make one punch on the opposite side, so you end up with 6 holes along the same edge.
- 5. Use scissors to trim bottom corners as shown in the Diagram 1-1.

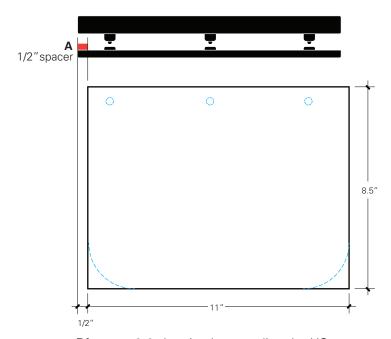


Diagram 1-1 showing how to align the US Letter clear plastic sheet with the hole punch to obtain the correct position for the first three holes, shown in blue.

See back for 6-hole shield template.

INSTRUCTION SHEET DESIGNED BY SYLVIA CEDEÑO - WWW.POGODESIGN.US

